BIOTECK®

CHondroGrid

FOR THE WELL-BEING OF MOVEMENT

Low molecular weight hydrolysed collagen



PREPARATION OF PRODUCT

PREPARATION OF CHONDROGRID® TO GUARANTEE STERILITY OF THE PRODUCT

WHAT IS NEEDED:

- A package of ChondroGrid®


DOUBLE STERILE PACKAGE

The double sterile package allows the product to be used even in the operating room




IN ADDITION, THE FOLLOWING ARE NEEDED:

- A syringe
- Two sterile needles
- A 2 ml bottle of sterile water

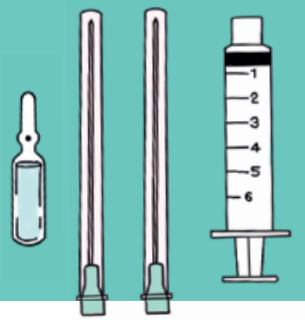


1

- Withdraw 2 ml of sterile, apyrogenic water into a clean 5 ml syringe fitted with a new needle.
- At the end of the procedure, protect the needle with its cap.

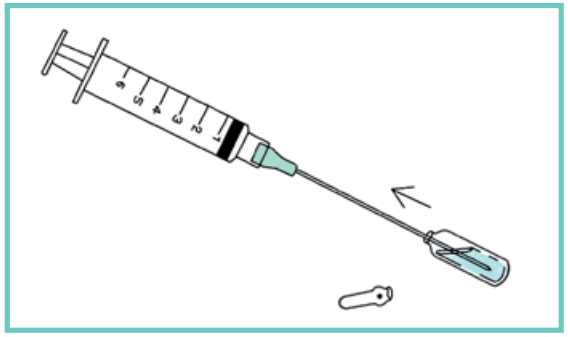


2

- Remove the ChondroGrid® bottle from its blister, avoiding touching the rubber stopper.
- Take the syringe previously loaded with 2 ml of water and inject the full contents into the ChondroGrid® bottle.
- Gently shake the bottle to help complete the ChondroGrid® suspension.

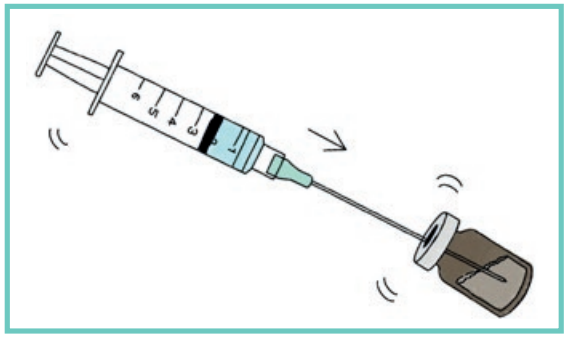

3


- Keeping the bottle tilted, completely aspirate the contents with the syringe.
- At the end of the procedure, protect the needle with its cap.

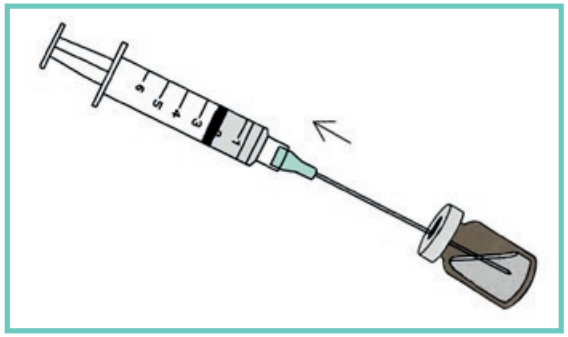

4


- Replace the first needle used for the preparation with the second needle to be used for injection.

PRESERVATION OF THE STERILITY CHAIN DURING PRODUCT HANDLING

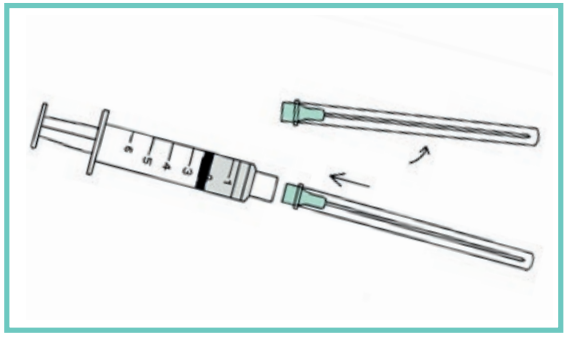

5

- After appropriately preparing the area, proceed with the injection of ChondroGrid®

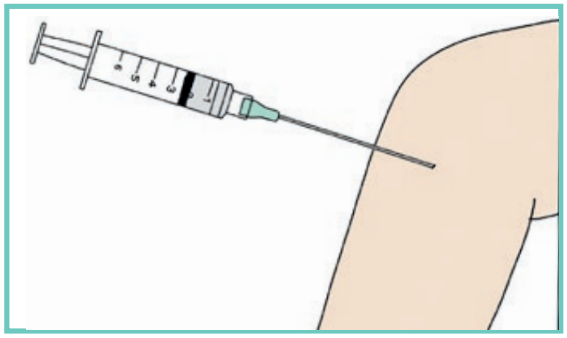

WATCH THE FULL VIDEO ON THE PREPARATION OF CHONDROGRID®



BIOTECK®.
INNOVATING BIOMATERIALS.

Bioteck® is an Italian company that has been producing bone substitutes and protective membranes a solution successfully used in orthopaedics, neurosurgery and oro-maxillo facial surgery since 1995.
Scientific research and innovation are the guiding principles that have enabled **Bioteck®** to patent new production processes and to create unique bio-materials of high quality in terms of performance level and safety guarantees. Materials now used in 72 countries worldwide.
At its multi-functional centre for research and development and thanks to state-of-the-art production processes, every day **Bioteck®** works to pursue its key objective: to innovate bio-materials.

WWW.BIOTECK.COM

BIOTECK ACADEMY.
SCIENTIFIC COMMUNITY FOR THE CULTURE OF THE CONSCIOUS CHOICE.

Bioteck Academy is the innovative scientific community which promotes the circulation and sharing of knowledge in the field of tissue regeneration applied to dentistry, maxillo-facial surgery, orthopaedics and neurosurgery.
Established as a hub for the clinical and scientific expertise focused on by **Bioteck®** spanning twenty years of research, today it is an entity open to all professionals who decide to join and share their own surgical experience.

WWW.BIOTECKACADEMY.COM

IMQ CERTIFIED

BIOTECK®

BIOTECK® Academy

Bioteck S.p.A.
Legal and Administrative Headquarters
Via E. Fermi, 49 – 36057
Arcugnano (VI) – Italy
Tel: +39 0444 289366
Fax: +39 0444 285272

Multifunctional Production, Research and Development Centre
Via G. Agnelli, 3 – 10020
Riva presso Chieri (TO) – Italy

© laziendinacreativa.it